

zensor)) Software

Two Part Software: zensoronline and zensor+




Data transmitted to server via
HTTP API



Device sends data to zensoronline via a simple HTTP API.
zensoronline is used to view the 30 sec event snapshots.

How it works

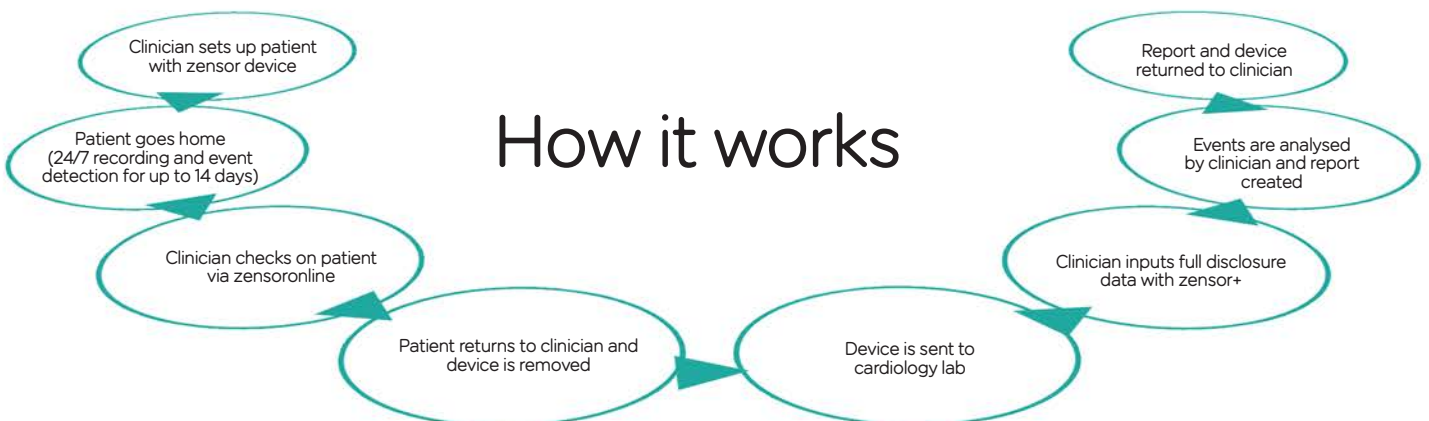


zensor+

Analysis Software

Device records full disclosure data during monitoring. The zensor+ software is then used to import this data for viewing by a clinician or cardiologist. All signals, beats, breaths and events are displayed and custom PDF reports created and saved.

zensor+ imports recorded full disclosure data from the device. Data is stored on the PC and user can view all beats, breaths, events and trend information

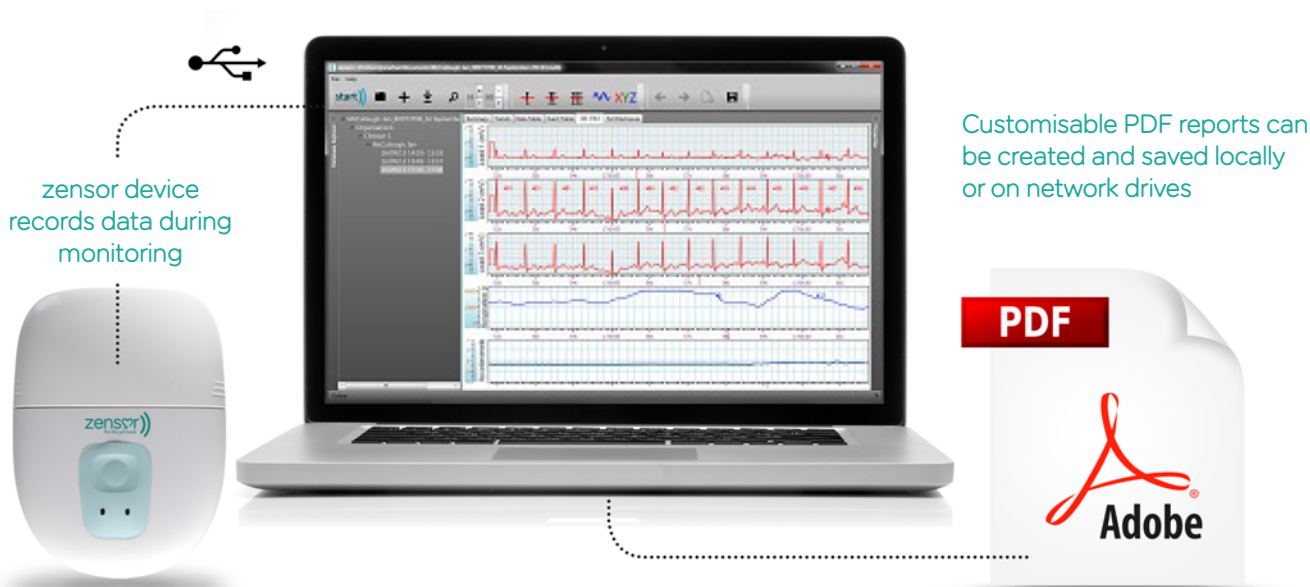

Feature	zensor+	zensoronline
Full Disclosure	✓	✗
Event Data	✓ (full disclosure pre and post event)	✓ (30s post-event)
Platform	Windows XP/7 (.NET 4)	Web based (PHP5/MySQL)
Near real-time monitoring	✗	✓ (events transmitted over network connection)*
Post monitoring	✓	✓ (event data remains available after monitoring)
CSV Export	Full disclosure	Event data
PDF Report	Customisable report creation	Daily report containing transmitted events

* Event data transmitted when network connection available

Let's Talk

Intelesens Limited
 17 Heron Road, Belfast,
 BT3 9LE,
 Northern Ireland.
 T: +44 (0)28 9073 6801
 F: +44 (0)28 9099 2433
 E: info@intelesens.com

www.intelesens.com